



IPL 5500
PURE SKIN PRO

beurer
beauty

IPL PURE SKIN PRO Long-lasting hair removal



Free "beurer MyIPL" app guides through the treatment and creates an individual treatment plan - for a professional IPL treatment at home



Up to 200,000 light pulses



With practical "auto-flash" mode for convenient use



Compact device for long-lasting hair removal - long-lasting results over several months

Compact design - perfect for treating small and hard-to-reach areas of the body

Around just 20 min. to treat the entire body¹⁾

3.1 cm² light surface for especially fast treatment

3 different intensity levels

Clinically tested technology for maximum safety at home

Maximum safety with 2-in-1 skin-type and skin-contact sensor

Up to 50% hair reduction after 3-4 treatments

State-of-the-art IPL light technology used by dermatologists for professional hair removal.

Deactivates the hair roots deep beneath the skin and in doing so decreases the chance of the hair growing back.

Integrated UV filter to protect the skin

Cannot be used on dark skin

May be ineffective on white or light-coloured hair



Mains operation:

Input: 100-240 V~; 50/60 Hz; 500 mA

Output: 12.0 V --- 1.8 A

Product measurements (LxWxH): 160 x 43 x 67 mm

Product weight: approx. 180 g

3 year guarantee

Sales unit: 4 / Shipping carton: -

EAN no.: 4211125 57514 4

Item no.: 575.14

* Dermatological study, 20 women, proderm GmbH, study no. 14.0309-11


¹⁾ at minimum power level for arms, legs, face, armpits and the bikini line




IPL 5500
PURE SKIN PRO

beurer
beauty


MODE OF ACTION

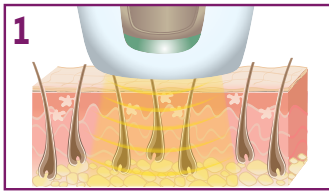

1
Melanin in the hair absorbs the light (IPL)

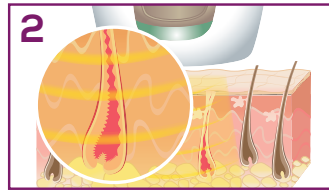


2
The warmth generated deactivates the hair follicle

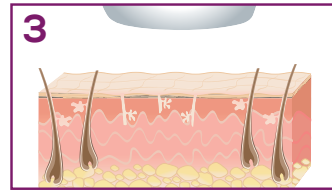

3
The hair falls out and does not continue to grow

HAIR AND SKIN COLOUR CHART

| | | Skin colour | | | | | |
|-------------|--|-------------|---|---|---|---|---|
| | | | | | | | |
| Hair colour | | X | X | X | X | X | X |
| | | ✓ | ✓ | ✓ | ✓ | ✓ | X |

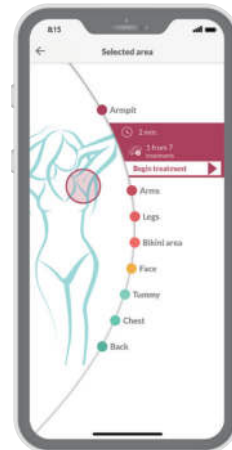

APPLICATION



works with

• iOS ≥ 12.0
• Android™ ≥ 8.0

List of supported devices

575.14_IPL5500_2021-09-13_DS_EN_BEU Subject to errors and changes

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.
Google Play and the Google Play logo are trademarks of Google LLC.
Android is a trademark of Google LLC.